


Classic T³™ Vaporizer

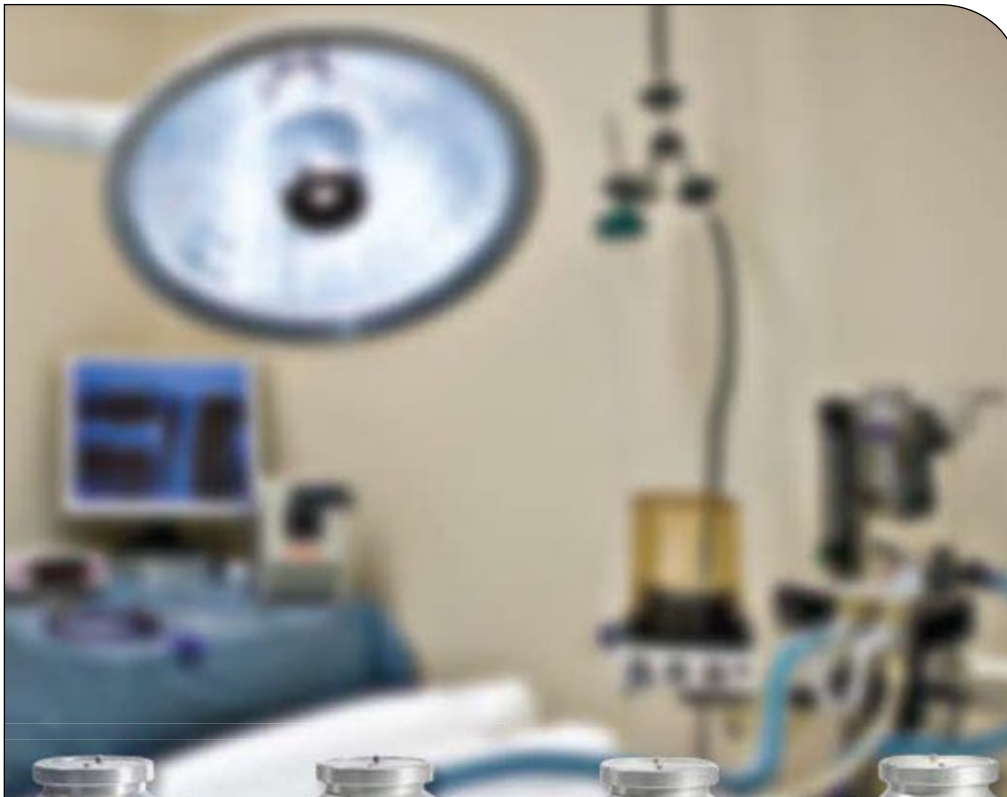

Classic T³™ Vaporizer

The Classic T³™ series is a precision style, out of circuit, single agent use vaporizer, that is temperature and flow compensated. Each vaporizer is clearly labeled for specific agent.

Advantages

- Isoflurane and Sevoflurane in both Funnel Fill and Key Fill
- All NEW parts, not refurbished
- Based on the original Tec 3 design
- One piece fill assembly with no solder joints, eliminating leak points
- Bold site glass for easy viewing of liquid volume
- Corrosion resistant
- RoHS compliant
- Three (3) year Warranty on Product and Calibration

VETERINARY

Classic T³ Vaporizer

Specifications

Inlet Flow Compensation

500 mls to 10 liters/min with minimal change in output concentration

Temperature Range Compensation

59° F to 95° F (15° C to 35° C)

The compensation assures relatively constant output with changes in ambient temperature and cooling due to vaporization

Calibration

69.8° F (21° C)

Calibration is warranted for three (3) years from date of shipment from the manufacturer

Accuracy

+/- 15% (within the above flow and temperature range specifications)

After three (3) year warranty period, annual calibrations are recommended

Duration of Use

Approximately 1.0 ml of liquid anesthetic is required to provide 200 ml of vapor

Liquid Capacity

125 ml liquid capacity
(35 mls is absorbed by a wick)

Services

- Exchange program available on all models
- Cleaning and calibration services available
- Verification service available during warranty period
- Contact customer service for additional information

Dimension and Weight

Funnel Fill

7.25" H x 5" W x 5.5" D
18.42 cm H x 12.7 cm W x 13.97 cm D
Weight: 16 lbs (7.27 kg)

Key Fill

7.25" H x 5" W x 6" D
18.42 cm H x 12.7 cm W x 15.24 cm D
Weight: 16 lbs (7.27 kg)

Catalog Number

Catalog #	Description
VCT302	Classic T ³ ™ Isoflurane Funnel Fill Vaporizer
VCT305	Classic T ³ ™ Sevoflurane Funnel Fill Vaporizer
VCT3K2	Classic T ³ ™ Isoflurane Key Fill Vaporizer
VCT3K5	Classic T ³ ™ Sevoflurane Key Fill Vaporizer



Easy viewing of liquid volume



One piece fill assembly

SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE
Please read the Instructions for Use supplied with the product for detailed instructions, warnings, and cautions.

Smiths Medical PM, Inc.

Veterinary

N7W22025 Johnson Drive, Waukesha, WI 53186 USA

Phone: 262-513-8500 Fax: 262-513-9069 Toll-Free USA: 888-745-6562

www.surgivet.com

Smiths Medical, part of the global technology business Smiths Group

Classic T³, SurgiVet and the Smiths Medical design mark are trademarks of the Smiths Medical family of companies. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners. © 2007 Smiths Medical family of companies. All rights reserved. Classic T³-Sell Sheet rev 01 10/07

smiths medical