

# SurgiVet® Equator® Convective Warming System



**The SurgiVet® Equator®** Convective Warming System is the leading choice in veterinary temperature management. The patented temperature control system ensures the desired therapy setting, while the high air flow promotes even distribution of warm air throughout the blanket and onto the patient...helping to improve patient outcomes and reduce the incidence of hypothermia.

For more information, visit [www.surgivet.com](http://www.surgivet.com).

**VETERINARY**

**SurgiVet**<sup>TM</sup>



Catalog #: WWVEQ5000RS

## SurgiVet® Equator® Convective Warming System

- Hose end temperature control automatically adjusts to deliver chosen air temperature—no matter how cold the environment
- Fast warm up and high air flow for quick and efficient thermal therapy
- Four temperature selections displayed on large digital readout for all therapy needs
- Over/under temperature and disconnect visual alarms on front panel ensures safety
- Compact size and lightweight (11.75" H x 9.5" W x 7.5" D, 15 lbs.)
- Can be mounted to an IV pole or set on the floor or tabletop
- Quiet air delivery
- Three year warranty-new
- One year warranty-remanufactured

### SurgiVet® Equator® Convective Warming System Includes:

- (1) 7 ft large air delivery hose
- (1) Full body blanket
- (1) Pediatric full body blanket
- (1) Lower body blanket
- (2) Pediatric underbody blankets
- Operator's manual

## Level 1® Snuggle Warm® Convective Warming Blankets

- Disposable and constructed from a soft fabric-like material
- Strong and durable—resists tears, punctures and fluids
- Designed to conform to the patient's body
- Provides even distribution of heat
- Large, heat-formed air holes in the blankets eliminate particulate material
- Latex free


Catalog #	Description	Blanket Size	Units/Case
WWVSW2001	Full Body Blanket	40" W x 80" L	10
WWVSW2002	Pediatric Full Body Blanket	40" W x 57.5" L	10
WWVSW2003	Upper Body Blanket	80" W x 40" L	10
WWVSW2004	Lower Body Blanket	40" W x 64" L	10
WWVSW2005	Neonate to Small Blanket (Requires small hose: 32491B13)	51" W x 28" L	10
WWVSW2006	Preemie to Neonate Blanket (Requires small hose: 32491B13)	40" W x 22" L	10
WWVSW2009	Small Pediatric Underbody Blanket	26" W x 37" L	10
WWVSW2011	Large Pediatric Underbody Blanket	40" W x 54" L	10



WWVSW2001




WWVSW2002

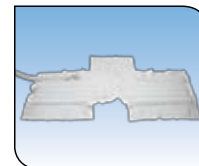

WWVSW2003



WWVSW2004

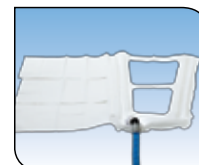

WWVSW2005



WWVSW2006



WWVSW2009



WWVSW2011

SPECIFICATIONS SUBJECT TO CHANGE WITHOUT NOTICE.

Please read the Instructions for Use supplied with the product for detailed instructions, warnings and cautions.

### Smiths Medical North America

N7W22025 Johnson Drive, Waukesha, Wisconsin 53186 USA  
 Phone: 262-513-8500 Fax: 262-542-0718 Toll-Free USA: 888-745-6562  
[www.surgivet.com](http://www.surgivet.com)

Smiths Medical, part of the global technology business Smiths Group

Equator, Level 1, Snuggle Warm, SurgiVet and the Smiths Medical design mark are trademarks of the Smiths Medical family of companies. The symbol © indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners. © 2009 Smiths Medical family of companies. All rights reserved. Equator® Rev.00 03/09

smiths medical